

Microfracture Awls for  
Chondral Defect Resurfacing



# Microfracture Awls

- 340mm length to access the deep hip or knee joint and facilitate leverage
- Wide range of tip angles (0, 30, 45, and 90 degrees)
- Ergonomic Design
- Create microfractures in cartilage and subchondral bone

# Microfracture Awls

The ConMed Linvatec ergonomically designed microfracture awl instrument set is designed to provide the surgeon with a full set of options for performing the hip and knee microfracture technique. The creation of localized trauma using this technique has been shown to enhance chondral resurfacing by providing a suitable environment for new tissue formation, taking advantage of the body's own healing potential.

Available in 10 different configurations, the awls are designed to work with commonly used cannulae, and are sufficiently long (340mm) to access the deep hip joint or knee joint and facilitate leverage. The awls are available in several different size and tip configurations.

## Light Awls


Designed to create microfractures in cartilage and debride tissue, available in 0, 30, 45 and 90 degree tip configurations. Curette also available.

## Heavy Awls

Designed to create microfractures in subchondral bone, available in 0, 30, 45 and 90 degree tip configurations. Features a strike plate on the proximal end and 3mm depth markings on the distal end (approximate depth needed to create local micro-fractures). Angled shaft heavy awl also available to facilitate access to difficult to reach areas.



Microfracture Instrumentation Set

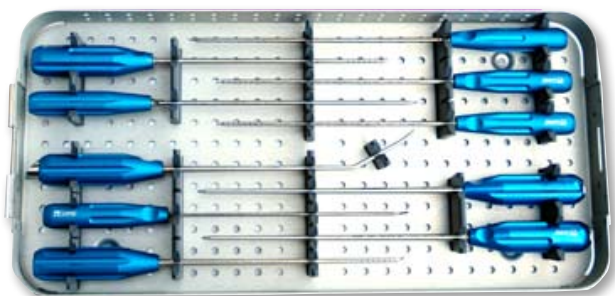

### Description

### Cat. No.

#### Light Awls

Microfracture Awl, 0 Degrees Light	MFX-0LT
Microfracture Awl, 30 Degrees Light	MFX-30LT
Microfracture Awl, 45 Degrees Light	MFX-45LT
Microfracture Awl, 90 Degrees Light	MFX-90LT
Microfracture Awl Curette, 45 Degrees	MFX-45CU

#### Heavy Awls

Microfracture Awl, 0 Degrees Heavy	MFX-0HV
Microfracture Awl, 30 Degrees Heavy	MFX-30HV
Microfracture Awl, 45 Degrees Heavy	MFX-45HV
Microfracture Awl, 90 Degrees Heavy	MFX-90HV
Microfracture Awl, Angled, 30 Degrees Heavy	MFX-30ANG

#### Sterilization

Microfracture Awl Sterilization Tray	MFX-TRAY
--------------------------------------	----------

### CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



11311 Concept Boulevard  
 Largo, FL 33773-4908  
 (727) 392-6464  
 Customer Service: 1-800-237-0169  
 FAX: (727) 399-5256  
 International FAX: +1 (727) 397-4540  
 email: customer\_service@linvatec.com  
 www.linvatec.com